

File D:\DATA\

Acq. Method : C:\CHEM32\2\METHODS\1-POS-3MIN.M

Analysis Method : C:\CHEM32\2\METHODS\1-POS-3MIN.M

Sample Info : Easy-Access Method: '1-POS-3min'

Method Info : Mobile Phase: A: water(10mM Ammonium hydrogen carbonate) B: ACN

Gradient: 5%-95% B in 1.6min


Flow Rate: 2.0ml/min

Column: XBridge C18,4.6*50mm,3.5um


Oven Temperature: 40 C

MS Range: 110-800


*ADC1 A, ELSD

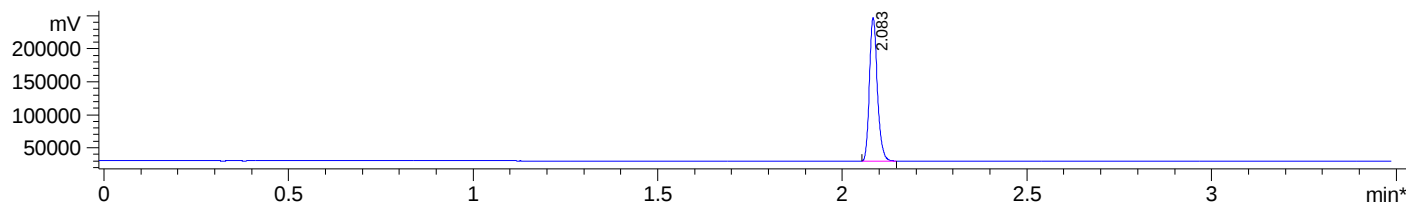


*DAD1 A, Sig=214,4 Ref=off

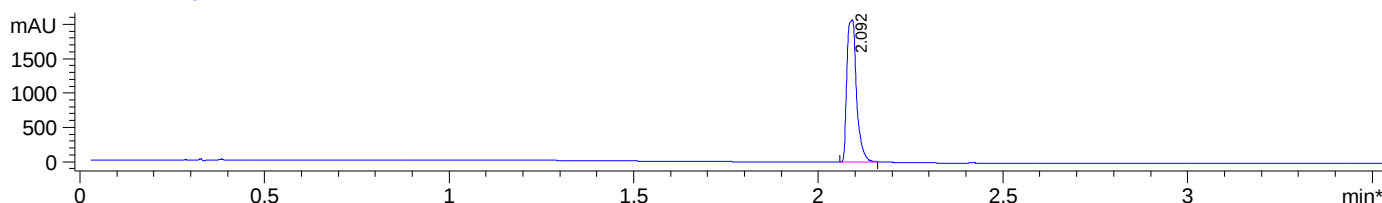

*DAD1 B, Sig=254,4 Ref=off

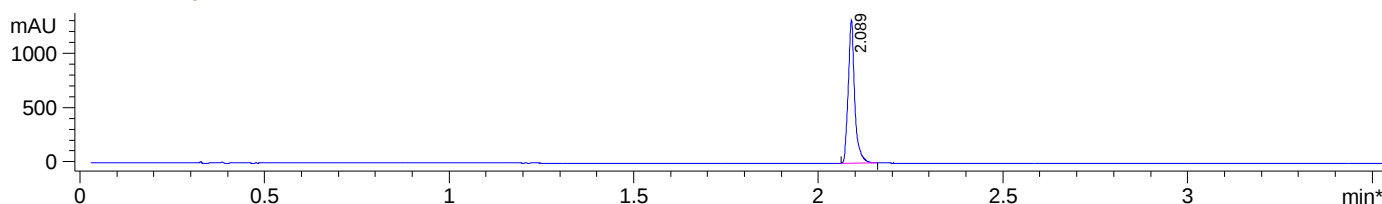

MSD1 TIC, MS File ES-API, Pos, Scan, Frag: 70

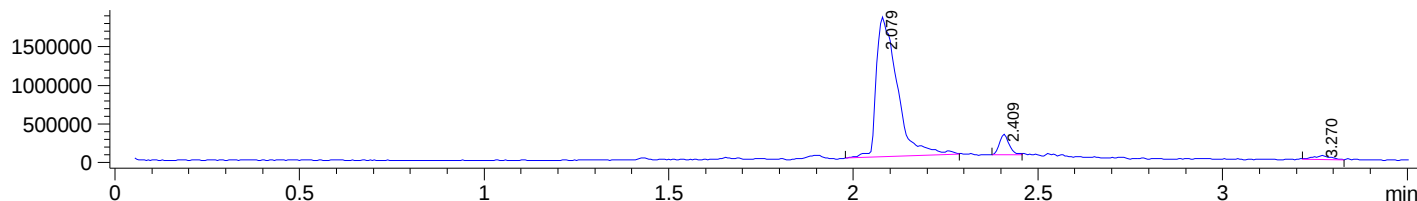

Integration Results for ADC1 A, ELSD

RetTim	Width	Area	Height	Area%
2.08	0.02	325996.72	217787.84	100.00

Integration Results for DAD1 A, Sig=214,4 Ref=off


RetTim	Width	Area	Height	Area%
2.09	0.02	3750.78	2072.61	100.00

Integration Results for DAD1 B, Sig=254,4 Ref=off


RetTim	Width	Area	Height	Area%
2.09	0.02	1625.73	1322.67	100.00

Ret. Time: 2.08

<<<< POSITIVE SPECTRA >>>>

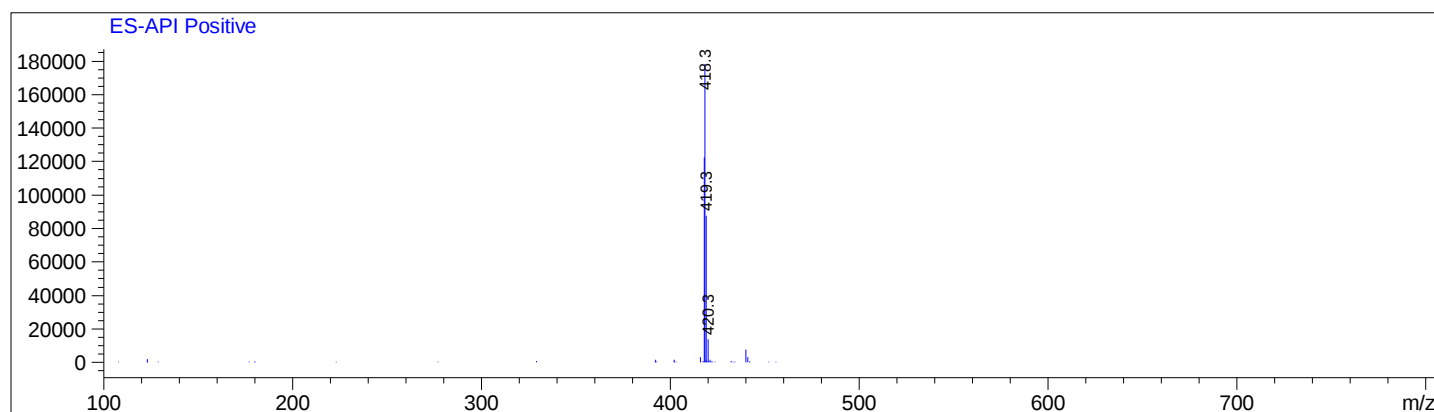

Ret. Time: 2.41

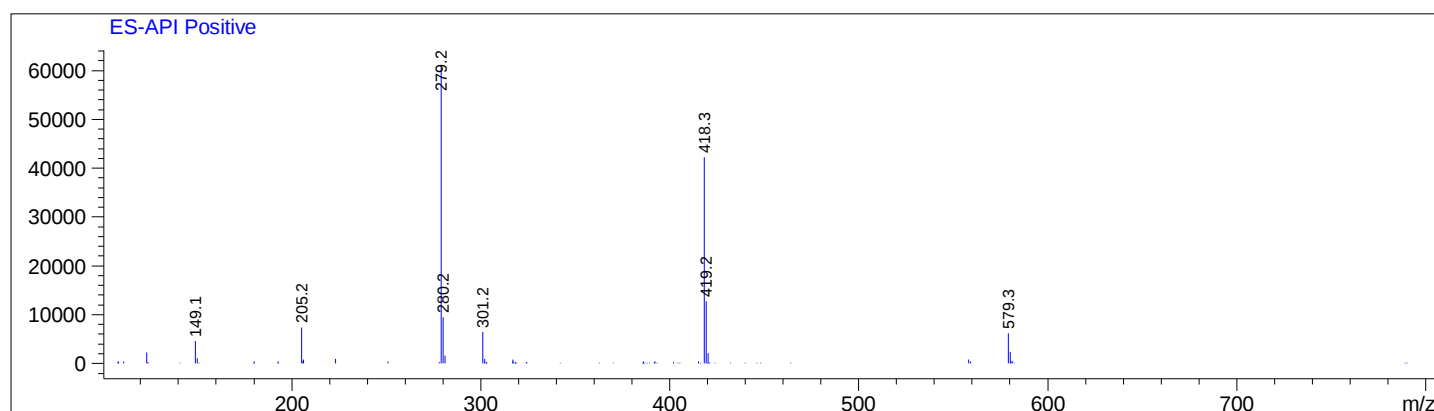

Ret. Time: 3.27

