

File D:\DATA\

Acq. Method : D:\METHODS\1-POS-3MIN.M

Analysis Method : D:\METHODS\1-POS-3MIN.M

Sample Info : Easy-Access Method: '1-POS-3MIN'

Method Info : Column:SunFire C18 (4.6x 50 mm,3.5 um)


Mobile phase: H2O(0.01%TFA) (A) / ACN(0.01%TFA) (B)

Elution program: Gradient from 5 to100% of B in 1.3 min at 2ml/min


Temperature: 45 °C

Detection: UV (214 , 4 nm) and MS (ESI, Pos mode ,110 to 1000 amu)


*ADC1 A, ELSD

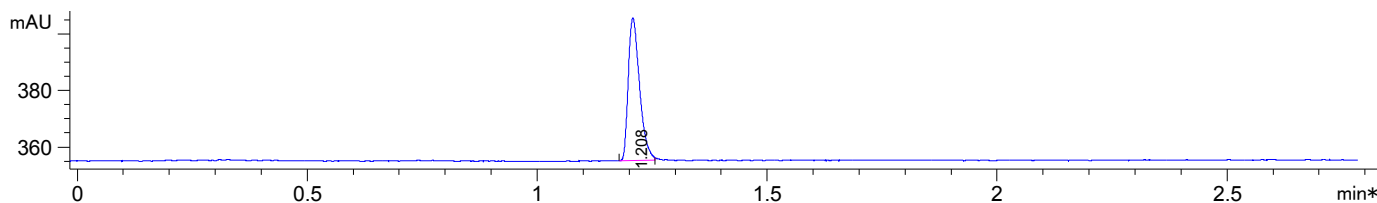


*DAD1 A, Sig=214,4 Ref=off

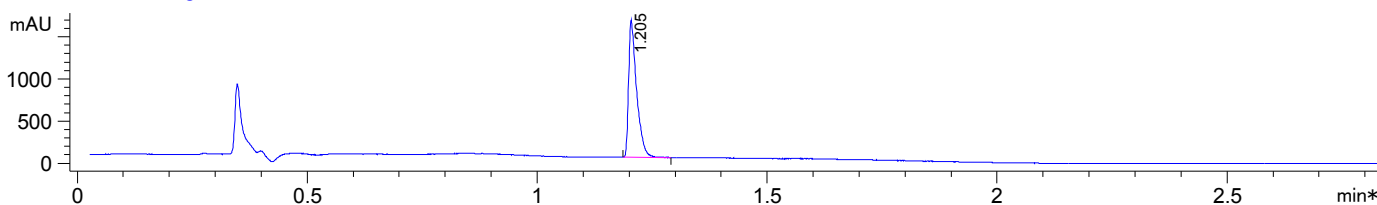

*DAD1 B, Sig=254,4 Ref=off

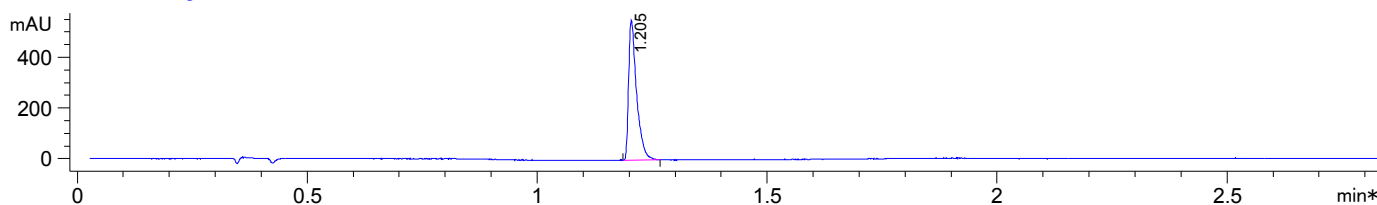

MSD1 TIC, MS File ES-API, Pos, Scan, Frag: 80

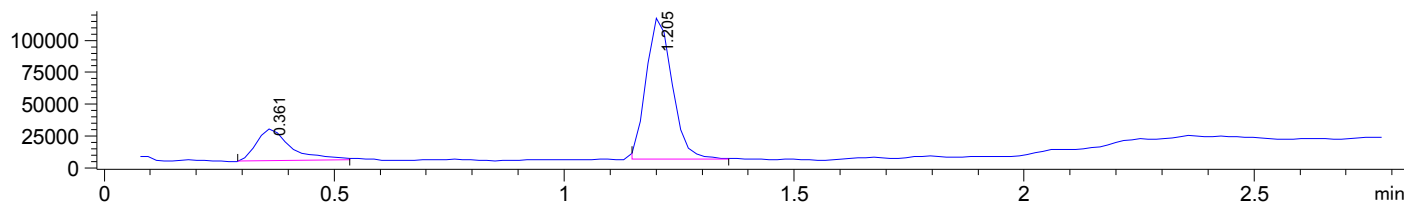

Integration Results for DAD1 A, Sig=214,4 Ref=off

RetTim	Width	Area	Height	Area%
1.20	0.02	1957.17	1604.27	100.00

Integration Results for DAD1 B, Sig=254,4 Ref=off

RetTim	Width	Area	Height	Area%
1.20	0.02	666.18	551.44	100.00

Ret. Time: 0.36 <<<< POSITIVE SPECTRA >>>> Ret. Time: 1.20

